

IUI Catheters

PRODUCT INFORMATION

Product no.	
MINICAT	1220905
FLEXCAT	1221805
GUIDECAT	1221905

MiniCat™

IUI-catheter

The small diameter of the Genitor reduces the residual volume, or dead space, to one tenth the volume of standard catheters. This is particularly desirable when performing inseminations with limited specimen volume. The small O.D. of this catheter also offers a high degree of flexibility and facilitates atraumatic access to the uterine fundus.

Residual volume: 0.05 ml.

FlexCat™

IUI-catheter with removable guidewire

Because of the inner “guide wire” (mandrel/stiletto) penetration of the cervix is easy – also if uterus is anteverted. The guide wire is removed before insemination. Although the larger O.D. of this catheter offers a small degree of rigidity to be of assistance in cannulating the cervix, its flexibility facilitates atraumatic access to the uterine fundus.

Residual volume: 0.38 ml


GuideCat™

IUI-catheter with fixed guidewire

This catheter includes both a separate tube for the sperm and a fixed “guide wire” in the catheter wall so that penetration of cervix is easy – also by anteverted uterus.

Residual volume: 0.35 ml.

Each catheter is wrapped in irradiation sterilized packs and delivered in packages of 25 catheters.

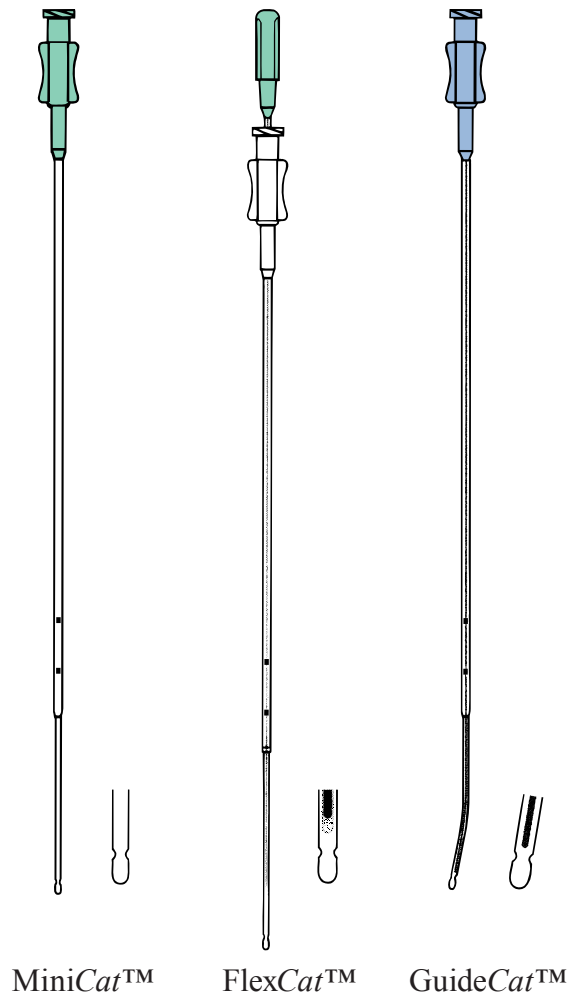

Vesterport 3 • 8000 Aarhus C • Denmark
Phone (+45) 86 76 06 99 • Fax (+45) 86 76 06 85
e-mail: info@cryosinternational.com
www.cryosinternational.com



Mouse
Embryo
toxicity
tested

MEA



Local distribution: